

# EDGE UTOPIA™

Enhanced Technology. Remarkably Beautiful. Ultra Secure Handle.

**Our Greatest Endodontic Systems Ever!**

Introducing Our New FireWire Blaze Heat Treatment

## EDGE X7 UTOPIA™

Glide  
Path

Initial  
Shaping

Optional  
Shaping



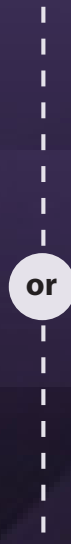
17/04



25/04



30/04



25/06



EDGEENDO®

Order Now

Free Shipping

EdgeUtopia.com

# EDGE UTOPIA™

## The Remarkable EdgeX7 Utopia™

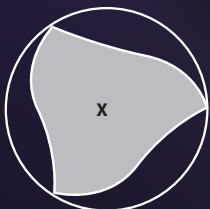
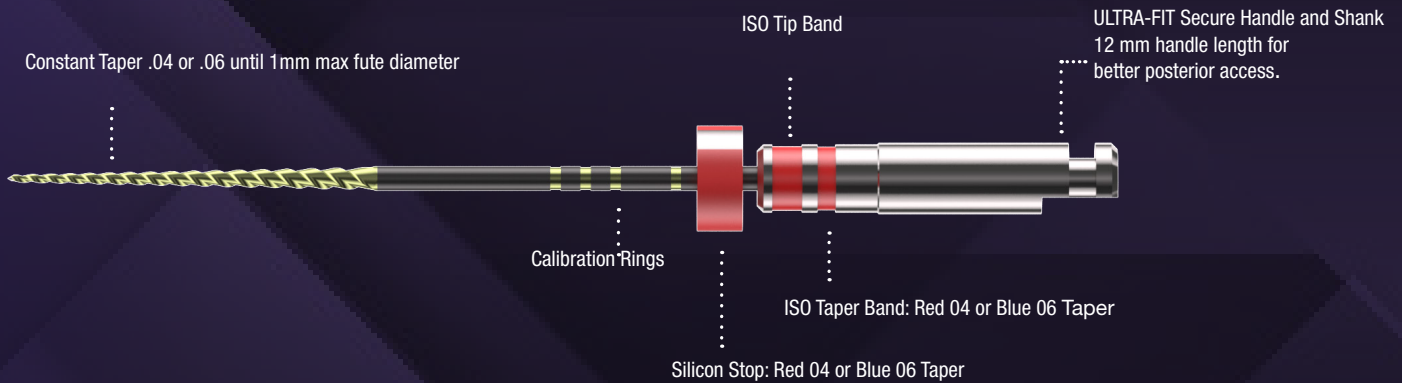
The New EdgeX7 Utopia™ delivers a Technology-Enhanced EdgeFile X7® with a Remarkably-Beautiful gold blade featuring our FireWire Blaze heat treatment. Plus, we added our new ultra secure, Ultra-Fit handle and shank.



### FireWire Blaze

FireWire Blaze—Years to develop and worth every minute. We wanted to outdo our original FireWire and with FireWire Blaze, we believe we have! We believe you will love the feel and flexibility as you negotiate the most challenging canals.

### Pre-Sterilized



Parabolic Cross Section

Available in 3 lengths: 21, 25 & 29 mm, 04 & 06 File Size

Maximum flute diameter 1 mm



Highly Flexible



No Bounce Back



High Strength



EDGEENDO®

Order Now

Free Shipping


EdgeUtopia.com

# EDGE X7 UTOPIA™

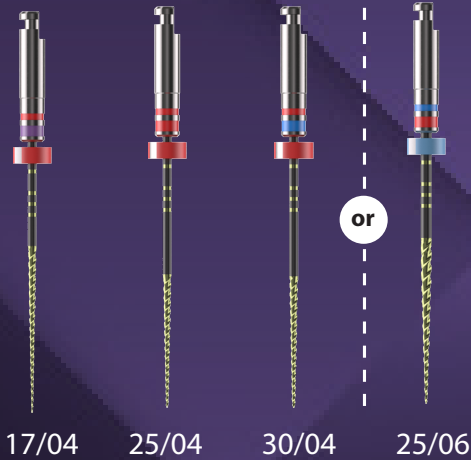Glide  
Path

Initial  
Shaping

Optional  
Shaping



Clockwise rotation  
Speed 300 to 500 rpm  
Torque setting 250 to 410 g/cm

Pre-Sterilized

## Simplified Technique Guide

1. Create straight line access to canal orifices.
2. Locate canals and explore using stainless steel hand instruments. Minimum size #10K file to working length recommended prior to rotary file use.
3. Irrigate before each hand or rotary file.
4. Use #17/04 file in one or more passes, alternating with small-sized hand files if necessary, until working length is reached. If more coronal flare is desired it can be achieved by incorporating the EdgeGlidePath file used in a brushing motion.
5. Next use #25/04 to working length passively; if instrument has not reached working length use additional shaping instrument #20/04 to working length.
6. If #25/04 reached working length with minimal resistance or if clinician desires a larger apical shape additional instruments can be used (#30, #35, #40, etc).

\*This technique is intended to be used as a guide only

**The New EdgeX7 Utopia™ is taking the original X7 to the next level of performance, cutting ability and flexibility. We are also excited to introduce our incredible FireWire Blaze NiTi along with our new ultra secure, Ultra-Fit handle and shank. This is our best X7 ever!**

**We know you will love the New EdgeX7 Utopia™...  
And so will your patients!**



EDGEENDO®

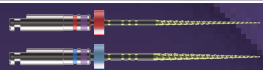

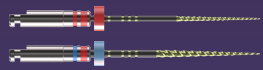
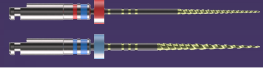
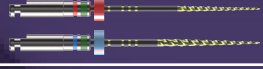


Order Now

Free Shipping

EdgeUtopia.com

# EDGE X7 UTOPIA™

- Available in .04 and .06 Constant Taper until 1mm max diameter
  - Variable Pitch
- Maximum flute diameter 1mm allows for minimally invasive preparation
- Parabolic Cross Section non cutting tip
  - Maximizes file cutting efficiency
- Reduced handle and shank length for increased posterior access
- ISO tip size 17-45
- Available lengths: 21, 25 & 29 mm
- FireWire Blaze technology
- Sterile Packaging
- Replaces EdgeFile® X7
- Same technique as EdgeFile® X7
- ULTRA-FIT Secure Handle and Shank
- Spectacular color coating finish
- Pre-Sterilized
- Free Shipping

| 6-Pack  | STERILE | Tip Size | Taper Size | Length        |               |             |
|---|---------|----------|------------|---------------|---------------|-------------|
|   |         |          |            | 21 mm         | 25 mm         | 29 mm       |
|  |         | 17       | 04         | X7U170421SE   | X7U170425SE   | X7U170429SE |
|   |         |          | 06         | X7U170621SE   | X7U170625SE   | X7U170629SE |
|  |         | 20       | 04         | X7U200421SE   | X7U200425SE   | X7U200429SE |
|   |         |          | 06         | X7U200621SE   | X7U200625SE   | X7U200629SE |
|  |         | 25       | 04         | X7U250421SE   | X7U250425SE   | X7U250429SE |
|   |         |          | 06         | X7U250621SE   | X7U250625SE   | X7U250629SE |
|  |         | 30       | 04         | X7U300421SE   | X7U300425SE   | X7U300429SE |
|   |         |          | 06         | X7U300621SE   | X7U300625SE   | X7U300629SE |
|  |         | 35       | 04         | X7U350421SE   | X7U350425SE   | X7U350429SE |
|   |         |          | 06         | X7U350621SE   | X7U350625SE   | X7U350629SE |
|  |         | 40       | 04         | X7U400421SE   | X7U400425SE   | X7U400429SE |
|   |         |          | 06         | X7U400621SE   | X7U400625SE   | X7U400629SE |
|  |         | 45       | 04         | X7U450421SE   | X7U450425SE   | X7U450429SE |
|   |         |          | 06         | X7U450621SE   | X7U450625SE   | X7U450629SE |
| Assorted Packs .04 (17, 20, 25, 30, 35)   |         |          |            | X7UASST0421SE | X7UASST0425SE | N/A         |
| Assorted Packs .06 (17, 20, 25, 30, 35)   |         |          |            | X7UASST0621SE | X7UASST0625SE | N/A         |

## GUTTA PERCHA POINTS

| Family               | Tip Size | Taper Size |
|----------------------|----------|------------|
| <b>EdgeFile® .04</b> |          |            |
| EdgeFile® GP .04     | 15       | 04         |
|                      | 20       | 04         |
|                      | 25       | 04         |
|                      | 30       | 04         |
|                      | 35       | 04         |
|                      | 40       | 04         |
|                      | 45       | 04         |
| <b>EdgeFile® .06</b> |          |            |
| EdgeFile® GP .06     | 15       | 06         |
|                      | 20       | 06         |
|                      | 25       | 06         |
|                      | 30       | 06         |
|                      | 35       | 06         |
|                      | 40       | 06         |
|                      | 45       | 06         |
| 50                   | 06       |            |



## PAPER POINTS

| Family               | Tip Size | Taper Size |
|----------------------|----------|------------|
| <b>EdgeFile® .04</b> |          |            |
| EdgeFile® PP .04     | 15       | 04         |
|                      | 20       | 04         |
|                      | 25       | 04         |
|                      | 30       | 04         |
|                      | 35       | 04         |
|                      | 40       | 04         |
|                      | 45       | 04         |
| 50                   | 04       |            |
| <b>EdgeFile® .06</b> |          |            |
| EdgeFile® PP .06     | 15       | 06         |
|                      | 20       | 06         |
|                      | 25       | 06         |
|                      | 30       | 06         |
|                      | 35       | 06         |
|                      | 40       | 06         |
|                      | 45       | 06         |
| 50                   | 06       |            |



**EDGEENDO®**

**Order Now**

**Free Shipping**

**EdgeUtopia.com**